TopClosure®

IVT Medical Ltd.

1SM RESPONDER KIT

Temporary wound closure for responders - ready when it matters.

WHAT'S INSIDE

TopClosure



STERILE

DRAPE STRIP



STERILE



Sterile Gauze Pad

10 x 10 cm



Elastic Bandage

7 cm x 2.7 m



2 Gloves - Size L



Alcohol Pad

70% Isopropyl Alcohol

STERILE



Alcohol Pad

70% Isopropyl Alcohol

STERILE

Single Use

Latex-Free

No Needles

Watch the application video

Scan before use

TopResponder.us

Application video

Temporary first-aid measure. Always seek medical attention as soon as possible.

See reverse side for instructions.

Six simple steps.



1

Clean


Wipe the skin around the wound with the antiseptic pads.



2

Cover

Place the sterile gauze pad over the wound.

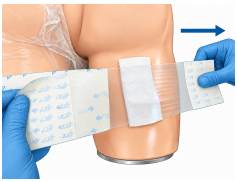

3

Seal


Apply the film drape over the gauze and surrounding skin. Smooth it down.



4

Attach

Stick the plates onto the drape, one on each side of the wound.

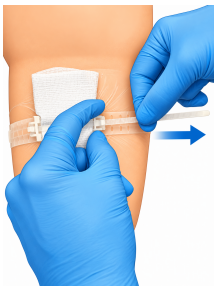

5

Close

Thread the strap and tighten gradually until the wound edges come together.



6

Secure

Wrap with the elastic bandage, keeping the strap end accessible.

Important safety information

The TopClosure® 1SM Responder Kit is a temporary first-aid measure for approximating and securing the edges of skin cuts and lacerations until professional medical care is available. It does not replace medical treatment. Seek medical attention as soon as possible after use.

Do not use for: deep or gaping wounds, heavy or spurting bleeding, puncture wounds, animal or human bites, wounds on the face or over joints, wounds showing signs of infection, or if you have a known sensitivity to adhesives or alcohol.

Single use only. Keep out of reach of children. Read all instructions before use.

TopClosure®

Application video

IVT Medical Ltd.